

0.035"



# FOXcross

PTA Balloon Dilatation Catheter

Totally coated  
balloon system

Low profile

Advanced shaft  
technology

Designed to get you there




## Jet coated balloon, shaft, tip and guide wire lumen for reduced friction

- Improved trackability
- Improved crossability, especially in calcified long lesions
- Reduced pull-back force

## Low profile – 5F sheath compatible up to 7 x 80 mm balloon

- Tapered tip
- Excellent rewrap
- High quality manufacturing

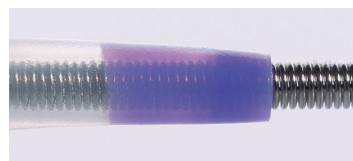


## Advanced shaft technology

- Crescent moon balloon lumen with increased cross section for rapid deflation
- Hydrophobic-coated guide wire lumen for excellent trackability

Crescent moon balloon lumen  
Increased cross section: **rapid deflation**

Hydrophobic-coated guide wire lumen:  
**excellent trackability**

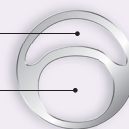

**Endovascular  
Commitment**

**Abbott**  
A Promise for Life



Designed to get you there

Recommended Guide Wire: 0.035"

## Product Information

80 cm Usable Shaft Length Ref. No.	135 cm Usable Shaft Length Ref. No.	Balloon Diameter [mm]	Balloon Length [mm]	Sheath Compatibility [F]	Nominal Pressure (NP) [atm]	Rated Burst Pressure (RBP) [atm]
10303-20	10312-20	3.0	20	5	10	18
10304-20	10313-20	4.0	20	5	10	18
10305-20	10314-20	5.0	20	5	9	17
10306-20	10315-20	6.0	20	5	8	16
10307-20	10316-20	7.0	20	5	8	16
10308-20	10317-20	8.0	20	6	9	15
10309-20	10318-20	9.0	20	6	9	14
10310-20	10319-20	10.0	20	7	7	13
10311-20	10320-20	12.0	20	7	5	9
10344-20	10345-20	14.0	20	7	5	8
10305-30	10314-30	5.0	30	5	9	17
10306-30	10315-30	6.0	30	5	8	16
10307-30	10316-30	7.0	30	5	8	16
10308-30	10317-30	8.0	30	6	9	15
10309-30	10318-30	9.0	30	6	9	13
10310-30	10319-30	10.0	30	7	7	13
10311-30	-	12.0	30	7	5	9
10303-40	10312-40	3.0	40	5	10	18
10304-40	10313-40	4.0	40	5	10	18
10305-40	10314-40	5.0	40	5	9	17
10306-40	10315-40	6.0	40	5	8	16
10307-40	10316-40	7.0	40	5	8	16
10308-40	10317-40	8.0	40	6	9	15
10309-40	10318-40	9.0	40	6	9	13
10310-40	10319-40	10.0	40	7	7	13
10311-40	10320-40	12.0	40	7	5	8
10344-40	10345-40	14.0	40	7	5	8
10304-60	10313-60	4.0	60	5	10	18
10305-60	10314-60	5.0	60	5	9	17
10306-60	10315-60	6.0	60	5	8	15
10307-60	10316-60	7.0	60	5	8	15
10308-60	-	8.0	60	6	9	15
10309-60	-	9.0	60	6	9	13
10310-60	-	10.0	60	7	7	13
10311-60	-	12.0	60	7	5	8
10303-80	10312-80	3.0	80	5	10	17
10304-80	10313-80	4.0	80	5	10	17
10305-80	10314-80	5.0	80	5	9	17
10306-80	10315-80	6.0	80	5	8	14
10307-80	-	7.0	80	5	8	13
10308-80	-	8.0	80	6	9	13
10304-100	10313-100	4.0	100	5	10	16
10305-100	10314-100	5.0	100	5	9	15
10306-100	10315-100	6.0	100	5	8	13
10307-100	10316-100	7.0	100	6	8	12
10304-120	10313-120	4.0	120	5	10	15
10305-120	10314-120	5.0	120	5	9	14
10306-120	10315-120	6.0	120	5	8	12
10307-120	10316-120	7.0	120	6	8	12

FOXcross is a trademark of the Abbott Group of Companies.

Contact your local representative for more information

### Abbott Vascular Australia:

Level 8, 15 Talavera Road, North Ryde, NSW 2113. Tel: 1800 550 939, Tel: +61 2 8879 2837, Fax: +61 2 8879 2805

### Abbott Vascular New Zealand:

Level 3, 6 Kingdom Street, New Market, Auckland. Tel: 0800 827 285, Tel: +64 9 524 6580, Fax: +64 9 524 6584

Photo(s) on file at Abbott Vascular. All drawings are artist's representations only and should not be considered as an engineering drawing or photograph. Products intended for use by or under the direction of a physician. Prior to use, it is important to read the package insert thoroughly for instructions for use, warnings and potential complications associated with use of this device. Information contained herein is for distribution for Australia and New Zealand ONLY.

For more information, visit our website at [abbottvascular.com](http://abbottvascular.com)

© 2011 Abbott. All rights reserved. ANZ 060069 10/11

